
Instructions for use

Dräger Alcotest 5000

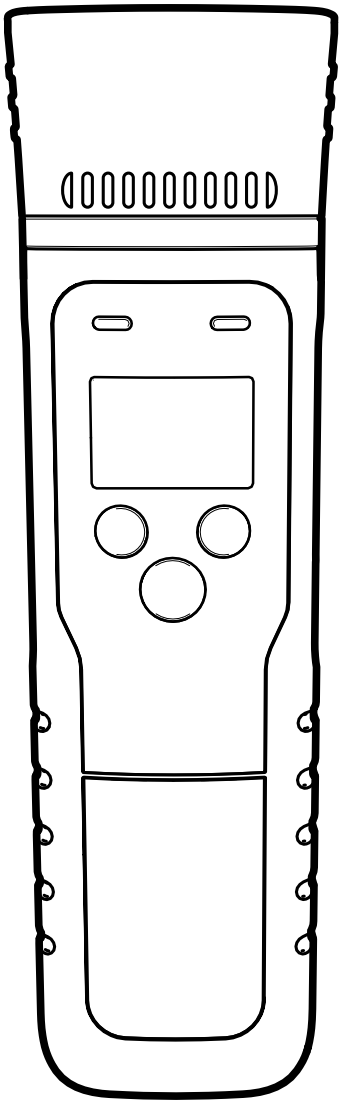More languages:
draeger.com/ifu

de, en, fr, es, ca, ptBR, it, nl, da, fi, no, sv, et, lv, lt, pl, ru, hr, sl, sk, sr, bs, cs, bg, uk, ro, hu, el, tr, ar, zh, zhTW, zhHK, ja, ko, hi, he

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1 Safety-related information 11

2 Conventions in this document 11

2.1 Meaning of the warnings 11
2.2 Typographical conventions 11

3 Description 11

3.1 Product overview 11
3.2 Intended use 11

4 Use 11

4.1 Prerequisites for use 11
4.2 Initial operation 12
4.3 Taking an automatic breath sample 12
4.4 Reading the result 12
4.5 Taking a manual breath sample (optional) 13
4.6 Taking a passive breath sample (optional) 13
4.7 Remove the Alcotest Hygienic-Flow funnel 13
4.8 Menu 13
4.9 Menu items 13

5 Troubleshooting 14

6 Maintenance 14

6.1 Replacing the batteries 14
6.2 Cleaning 14

7 Storage 14

8 Disposal 14

9 Technical data 15

10 Order list 15

1 Safety-related information

- Before using this product, carefully read these instructions for use and those of the associated products.
- Strictly follow the instructions for use. The user must fully understand and strictly observe the instructions. Use the product only for the purposes specified in the intended use section of this document.
- Do not dispose of the instructions for use. Ensure that they are stored and used appropriately by the user.
- Only trained and competent users are permitted to use this product.
- Do not use a faulty or incomplete product. Do not modify the product.
- Notify Dräger in the event of any component fault or failure.
- Comply with all local and national rules and regulations associated with this product.
- Only trained and competent personnel are permitted to inspect, repair and service the product. Dräger recommend a Dräger service contract for all maintenance activities and that all repairs are carried out by Dräger.
- Use only genuine Dräger spare parts and accessories. Otherwise, the proper functioning of the product may be impaired.

ℹ These instructions for use can be downloaded in other languages in electronic format from the appropriate product page (www.draeger.com/ifu) or obtained free of charge as a printed copy from Dräger or the specialist shop.

- Only use chargers, power packs, batteries, or rechargeable batteries that Dräger has approved for this product.

2 Conventions in this document

2.1 Meaning of the warnings

The following warnings are used in this document to alert the user to possible hazards. A definition of the meaning of each warning is as follows:

Warning symbols	Signal word	Classification of warning
⚠	WARNING	Indicates a potentially hazardous situation. If not avoided, it could result in death or serious injury.
⚠	CAUTION	Indicates a potentially hazardous situation. If not avoided, it could result in physical injury. It may also be used to alert against unsafe practices.
	NOTICE	Indicates a potentially hazardous situation. If not avoided, it could result in damage to the product or environment.

2.2 Typographical conventions

Text Text in bold identifies labels on the device and text on the display.

► This triangle in warnings identifies the options available for avoiding the hazard.

> The greater-than symbol identifies a navigation path in a menu.

ℹ This symbol indicates information that may make the use of this product easier.

3 Description

3.1 Product overview

See fold-out page:

1	Alcotest Hygienic-Flow funnel
2	Alcotest Hygienic-Flow funnel holder
3	LEDs (red, green)
4	Display
5	button (open menu; set value)
6	button (switch instrument on/off; confirm entries)
7	Wrist strap
8	button
9	USB port (not for charging the device, for service personnel only)
10	Eye for wrist strap

3.2 Intended use

The Dräger Alcotest® breath alcohol testing device[1] 5000 is used to measure the concentration of alcohol in human breath, in the air above substances or in the ambient air (passive test). This measurement is used to determine a person's blood alcohol level. The detection results should not be used as grounds for medical or therapeutic treatment or for legal actions.

4 Use

4.1 Prerequisites for use

NOTICE
Avoid sampling in a draught. It could falsify the test result.

- After the last alcohol consumption, wait for approx. 15 minutes.
- The wait time cannot be shortened by rinsing the mouth with water or a non-alcoholic drink.
- Residual mouth alcohol or aromatic drinks (e.g. fruit juice), alcoholic mouth sprays and medicines can falsify the measurement.
- Burping and vomiting can falsify the measurement result.
- A draught can falsify the measurement result.

4.1.1 Before sampling

The test person must breathe regularly and calmly before sampling. The test person must be capable of providing the minimum breath volume required for sampling.

1) Alcotest® is a registered trademark of Dräger.

4.1.2 Attaching the wrist strap

Use the wrist strap to secure the device to your wrist.

1. Feed the fastening lug of the wrist strap through the eye on the device.
2. Pull the top part of the wrist strap through the fastening lug.

4.1.3 Attaching telescopic rod

To be able to reach test persons who are out of reach (such as a driver sitting in the cab of his HGV), the device can be connected to a telescopic rod using a 1/4″ connecting piece.

- A telescopic probe is not included in the scope of delivery of the device.
- Adhere to the instructions for use provided by the manufacturer of the telescopic rod.

4.1.4 Preparing for use

1. Make sure that the batteries are charged and correctly inserted.
2. Make sure that a clean Alcotest Hygienic-Flow funnel has been fitted.
3. Make sure that the instrument is securely and safely held by correct use of the wrist strap or the telescopic rod.

4.1.5 Turning on the device

1. Press and hold the OK button for approx. 1 second.
2. Wait until the display indicates that the device is ready to take a measurement.

4.1.6 Turning off the device

- Press and hold the OK button for approx. 3 seconds.

ℹ The device will turn itself off 4 minutes after last being operated.

4.2 Initial operation

When the device is turned on for the first time, a setup wizard appears. No setup wizard appears for pre-configured devices.

The following settings are made using the setup wizard:

- Language
- Date
- Time format
- Time

All settings can be subsequently changed via the menu.

ℹ Loss of default settings
If the battery is removed for longer than 5 minutes, then the date and time will need to be reset.

4.3 Taking an automatic breath sample


⚠ WARNING
Hygiene
Distorted measurement results

- ► Do not touch the top part of the Alcotest Hygienic-Flow funnel.
- ► Always maintain a distance of 2–3 cm from the subject or object to be tested.
- ► Do not touch the skin of a test person and secretions or surfaces with the Alcotest Hygienic-Flow funnel.
- ► Avoid direct contact with substances to be tested.
- ► Replace the contaminated Alcotest Hygienic-Flow funnel.

⚠ WARNING
Damage to the device
If several tests are carried out in a short period of time, condensation can collect inside the Alcotest Hygienic-Flow funnel which can then get into the device.

- ► Replace the Alcotest Hygienic-Flow funnel in good time.

1. Remove the Alcotest Hygienic-Flow funnel from its packaging. Touch the Alcotest Hygienic-Flow funnel using the gripping areas on the sides and mount it firmly on the top of the device (only possible in one direction).
2. Ask the subject to blow evenly and without interruption into the Alcotest Hygienic-Flow funnel without touching it. An adequate breath flow is signalled by a continuous tone sounding.
 While the sample is being taken, the display shows a progress bar (see images **A** and **B**).
 Once the sampling has been successful, the tone ends with an audible click.
 The display shows a spinning progress bar during an ongoing analysis process.
 The result is shown on the display after a few seconds and is stored in the instrument.

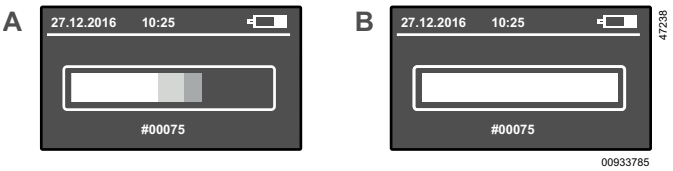4.4 Reading the result

After a few seconds (depending on the concentration measured and the temperature of the device), the test results are shown on the display and saved to the instrument. The following 2 different results are possible:

C) No alcohol has been found to be present. The test is passed. In addition, the green LED lights up. After a few seconds, the instrument is ready to use again.

D) Alcohol has been found to be present. The test is not passed. In addition, the red LED lights up.

1. Press the OK button to continue with the next measurement.

- During the sensor recovery time, the display shows a rotating progress symbol.
- Repeat the process (see "Taking an automatic breath sample")

Confirmation of the test results ***Alcohol*** and ***No alcohol*** can be configured in the menu as follows:

- The test results can be confirmed by pressing the OK button.
 Or
- The test results are confirmed automatically after 2 seconds.

4.5 Taking a manual breath sample (optional)

If the test person is not capable of providing the minimum volume required for sampling, sampling can be triggered manually.

1. Follow steps 1 to 3 of the automatic breath sample.
2. If the progress bar is shown after the breath sample has started, briefly press the OK button.
 When the results are displayed, they also show ***MANUAL***.

4.6 Taking a passive breath sample (optional)

With passive sampling, sampling is triggered manually by the person using the device. Possible applications for this include testing the ambient air or a person's end-expiratory breath for the presence of alcohol.

1. Press and hold the OK button for approx. 1 second until the startup screen appears.
2. Hold the prepared device **without** the Alcotest Hygienic-Flow funnel in the air to be analysed.
3. Briefly press the OK button to trigger the measurement.
 When the results are displayed, they also show ***PASSIVE***.

4.7 Remove the Alcotest Hygienic-Flow funnel

- Hold the Alcotest Hygienic-Flow funnel using the grip areas on its sides, pull it out of the device in an upwards motion and dispose of it.

4.8 Menu

To access the menu:

- Press the ▼ button.

To navigate to the menu or to configure the settings:

1. Press the ▼ or ▲ button.
2. Press the OK button to confirm the feature/input chosen.

To exit the menu:

- Choose the menu item **Back** in the sub-menus.
- Choose **Analyse** in the main menu.

4.9 Menu items

The menu contains the following items:

4.9.1 Measurement results

Shows all saved results with measurement unit, date and time, as well as sample number and the information as to whether the sample was taken manually or passively.

4.9.2 Infoscreen

Shows the part number and the serial number of the device, as well as the serial number of the EC sensor.

4.9.3 Next service

Displays the time until adjustment or service is due in days.

4.9.4 Date and time

Sets the date and time.

1. Change the digits using the ▼ or ▲ button.
2. Confirm with the OK button.

4.9.5 Contrast

Adjusts the contrast.

1. Change the contrast setting using the ▼ or ▲ button.
2. Confirm with the OK button.

4.9.6 Language

Sets the language.

1. Select the language using the ▼ or ▲ button.
2. Confirm with the OK button.

4.9.7 Software version

Displays the version number and part number of the software.

4.9.8 Temporary counter

Resets the test counter.

1. Reset the counter by confirming with the OK button.

4.9.9 Display

Adjusts the brightness of the display screen.

1. Set the brightness using the ▼ or ▲ button.
2. Confirm with the OK button.

4.9.10 Loudspeaker

Turns the volume on or off.

1. Select the desired feature (on or off) using the ▼ or ▲ button.
2. Confirm with the OK button.

4.9.11 Automatic shut-down

Sets the time until automatic shut-down.

1. Set the time using the ▼ or ▲ button.
2. Confirm with the OK button.

4.9.12 Signing in

For service personnel only.

1. Enter the four-digit code to open the service menu.

4.9.13 Automatically to the next test

Sets whether the instrument should automatically change to the next test after breath sampling. The setting can be set for both measurement results.

1. Select the (***Alcohol*** or ***No alcohol***) measurement result with the ↓ or ↑ button.
2. Confirm with the OK button.
3. Select with the ↓ or ↑ or ***activated*** or ***deactivated*** button.
4. Confirm with the OK button.

4.9.14 Turning off

Turns off the device.

1. Briefly press the OK button.

✓ The device turns off.

5 Troubleshooting

In the event of a device malfunction, an error message will be shown. The number which appears underneath the message is used for service functions. If the malfunction does not disappear after turning the device on and off several times, contact DrägerService.

Malfunction/Error	Cause	Remedy
Device cannot be turned on.	The batteries are empty.	Change the batteries.
Volume too small/sample invalid	The test person is blowing unevenly or not strongly enough.	Blow into the Alcotest Hygienic-Flow funnel with more force and without interruption.
Device outside temperature range	The device is too hot or too cold.	Allow the device to cool down or warm up.
↑↑↑ (measuring range exceeded)	The measuring range has been exceeded. Residual mouth alcohol can lead to incorrect measured values.	At least 15 minutes wait time must have passed since alcohol was last consumed through the mouth. To repeat, press the OK- button.
Sample interruption	Irregular exhalation, abrupt ending or inhalation	Repeat breath sample by pressing the OK button
Error number	Device error	Turn the device off and then back on. If the error persists, contact DrägerService or your retailer.

6 Maintenance

The device is delivered with default adjustment settings. The date of the last adjustment can be recalled in the menu.

Dräger recommends that the device is adjusted every 24 months. Adjustment must only be carried out by authorized persons, e.g. DrägerService.

Calibration certificates:
https://www.draeger.com/productioncertificates

ℹ Observe the national guidelines for adjustment.

6.1 Replacing the batteries

1. Turn off the device.
2. Push off the battery cover.
3. Change the batteries (3 x 1.5 V AA alkaline batteries). Ensure that the polarity is correct when doing so.
4. Reattach the battery cover.
5. Check whether the time has been set correctly.

ℹ If changing the battery takes more than 5 minutes, the time needs to be reset. Measurements cannot be performed if the time has not been set.

6.2 Cleaning

Clean the device using an ethanol-based cleaning agent.

⚠ WARNING
Explosion hazard
Alcohol can ignite.
► Do not use cleaning agents containing alcohol near potential ignition sources.

To clean the device:

1. Apply the cleaning agent to a disposable cloth and wipe the device. Do not apply cleaning agent directly to the device.

ℹ Make sure that no liquid enters the gas inlet.

2. Carry out a passive sampling to ensure there is no residual ethanol inside the device and that the ***No alcohol*** result is displayed.
3. If ***Alcohol*** continues to be shown, repeat the passive sampling after a short break (10-15 minutes) until***No alcohol*** is displayed.
4. If ***Alcohol*** continues to be shown, store the device for up to 24 hours to allow the cleaning agent to evaporate. Repeat the passive sampling after some time.

Information about the disinfection of Alcotest devices: Document 9300698 at www.draeger.com/IFU.

7 Storage

⚠ WARNING
Damage to the device
Leaking battery acid can cause damage to the device.
► Remove the batteries if the device is not used for a prolonged period of time.

8 Disposal

This product must not be disposed of as household waste. This is indicated by the adjacent symbol.
You can return this product to Dräger free of charge. For information please contact the national marketing organizations or Dräger.

Batteries must not be disposed of as household waste. They are therefore marked with the adjacent symbol. Dispose of batteries at battery collection centres as specified by the applicable regulations.

9 Technical data

Measuring principle	
Electrochemical sensor	
Ambient conditions	
During operation	-15 °C to +50 °C (+5 °F to 122 °F) 10 to 100 % RH (non-condensing) 600 to 1300 hPa
During storage	-20 °C to +60 °C (-4 °F to +140 °F) 15 to 75 % RH (10 to 100 % RH <48 h)
Measurement results	
0 mg/L to 0.029 mg/L	no alcohol (✓) displayed
0.030 mg/L to 2.5 mg/L	alcohol (✗) displayed
More than 2.51 mg/L	Over range (↑↑↑) displayed
Sensitivity drift	
Typically 0.4 % of the measured value per month	
Wait time for measurement readiness following a previous measurement of:	
0.00 mg/L, 2 s	0.00 ‰, BAC, 2 s
0.25 mg/L, 20 s	0.50 ‰, BAC, 20 s
0.5 mg/L, 40 s	1.00 ‰, BAC, 40 s
1.0 mg/L, 80 s	2.00 ‰, BAC, 80 s
>1.50 mg/L, 120 s	>3.00 ‰, BAC, 120 s
Adjustment	
Recommendation: every 24 months	
Dimensions	
63 mm x 219 mm x 41 mm (W x H x D)	
Weight	
approx. 245 g	
Power supply	
3 x 1.5 V AA alkaline batteries	
The batteries are designed for up to 5,000 measurements (depending on use and ambient conditions)	
Backup	
The last 500 results are saved.	
Interface	
Micro-USB port, for service personnel only	
CE marking	
Electromagnetic compatibility	

10 Order list

Accessories	Part no.
Pack with Alcotest Hygienic-Flow funnels, 10 pcs.	82 27 718
Alkaline battery, Mignon, LR6, AA, 1.5 V	13 35 804
Reflective tape set, yellow, Alcotest® 5000	83 27 742
Reflective tape set, orange, Alcotest® 5000	83 27 743
Calibration adapter, Alcotest® 5000	83 27 741

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© Dräger Safety AG & Co. KGaA
Dräger Safety AG & Co. KGaA
Draeger Safety UK Ltd