
Instructions for use

Dräger Alcotest® 7000

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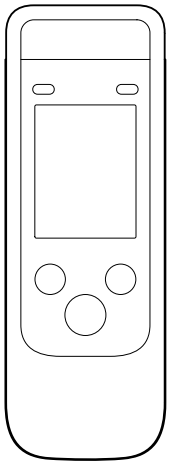Contents

1 Safety information 32

2 Conventions in this document 33
2.1 Meaning of the warning notes. 33
2.2 Typographical conventions 33

3 Description 33
3.1 Product overview 33
3.2 Intended use 34

4 Use 34
4.1 Starting the device for the first time 34
4.2 Menu and context menu 35
4.3 Prerequisites for Use 36
4.4 Preparations for use 36
4.5 Changing the measuring mode. 37
4.6 Breath sample. 37
4.7 Breath sampling with the funnel 40
4.8 Turning off the device 41
4.9 Removing the mouthpiece 41
4.10 Removing the funnel. 41
4.11 Menu items 42
4.12 GPS (Global Positioning System, optional) 45
4.13 Printing Results. 45
4.14 Exporting data. 46
4.15 App connection 46

5 Troubleshooting 47

6 Maintenance 48
6.1 Charging the rechargeable battery 49
6.2 Replacing the Batteries. 49
6.3 Cleaning 50

7 **Storage** 50

8 **Disposal** 51

9 **Technical data** 51

10 **Order list** 54

1 Safety information

- Before using this product, carefully read these instructions for use and those of the associated products.
- Strictly follow the instructions for use. The user must fully understand and strictly observe the instructions. Use the product only for the purposes specified in the intended use section of this document.
- Do not dispose of the instructions for use. Ensure that they are stored and used appropriately by the user.
- Only trained and competent users are permitted to use this product.
- Do not use a faulty or incomplete product. Do not modify the product.
- Notify Dräger in the event of any component fault or failure.
- Comply with all local and national rules and regulations associated with this product.
- Only trained and competent personnel are permitted to inspect, repair and service the product. Dräger recommends a Dräger service contract for all maintenance activities and that all repairs are carried out by Dräger.
- Use only genuine Dräger spare parts and accessories. Otherwise, the proper functioning of the product may be impaired.
- Only use battery chargers, power packs, batteries, or rechargeable batteries that Dräger has approved for this product.
- The power supply unit 3701090 has not been approved for maritime application.

These instructions for use are available in further languages and can be downloaded in electronic form on the according product website (www.draeger.com) or a printed version can be ordered free of charge via Dräger or specialist shops.

2 Conventions in this document

2.1 Meaning of the warning notes

The following warning notes are used in this document to notify users of possible dangers. The meanings of the warning notes are defined as follows:

Alert icon	Signal word	Classification of warning
⚠	WARNING	Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.
⚠	CAUTION	Indicates a potentially hazardous situation which, if not avoided, could result in injury. It may also be used to alert against unsafe practices.
	NOTICE	Indicates a potentially hazardous situation which, if not avoided, could result in damage to the product or environment.

2.2 Typographical conventions

Text Text marked in bold denotes labeling on the device and on-screen messages.

► This triangle labels possible methods of avoiding the hazards mentioned in warning notices.

> The greater-than sign denotes the navigational path in a menu.

This symbol denotes information that make using the product easier.

3 Description

3.1 Product overview

See Figure 1 on the fold-out page:

1 Breath alcohol testing device

2	USB-C port
3	Wrist strap
4	button (set value, open context menu)
5	Display
6	Mouthpiece
7	Brightness sensor
8	LEDs
9	button (load menu; set value)
10	button (turn the device on or off; confirm input)
11	Mouthpiece receptacle/funnel receptacle
12	Funnel

3.2 Intended use

The Dräger Alcotest® 7000 breath alcohol testing device[1)] is used to measure the concentration of alcohol vapor in a person's breath. This measurement is then used to determine a person's blood alcohol level. The measurement should not be used as grounds for medical or therapeutic treatment.

4 Use

4.1 Starting the device for the first time

4.1.1 Attaching the wrist strap

Use the wrist strap (3) to secure the device to the wrist (see Figure 1 on the fold-out page).

1. Thread the fastening lug of the wrist strap through the hole in the casing and open up the loop.
2. Pull the lower part of the wrist strap through the loop that has just been made.

1) Alcotest® is a registered trademark of Dräger.

4.1.2 Charging the rechargeable battery

Observe the following information: "Charging the rechargeable battery", page 49.

4.1.3 Inserting batteries (for battery-operated devices)

Observe the following information: "Replacing the Batteries", page 49.

4.1.4 Turning on the device/setup wizard

The brightness of the display automatically adapts to ambient conditions.
If the brightness sensor (7) - see Figure 1 on the fold-out page - is covered over, the device dims the brightness of the display, making it difficult to read the display in bright light conditions.
Never cover the brightness sensor with labels etc.

- Press the OK button for approximately 2 seconds to turn on the device.

⇒ The setup wizard opens.
This step does not apply in the case of preconfigured devices.

The setup wizard will configure the following settings:

- Language
- Date format
- Date
- Time format
- Time
- Measurement module

All settings can be altered at any time via the menu.

4.2 Menu and context menu

To access the menu:

- Press the button.

To access the context menu:

- Press the button

To navigate in the menu or context menu or to configure the settings:

1. Press the or button.
2. Press the OK button to confirm the selected feature/entry.

To exit the menu or context menu:

- Select the ***Back*** menu entry.

4.3 Prerequisites for Use

The device is adjusted in the factory before delivery. Dräger recommends that the device be adjusted every 12 months by DrägerService.

- Wait at least 15 minutes after consuming alcohol.
- This wait time cannot be shortened by rinsing the mouth with water or non-alcoholic drinks.
- Residual mouth alcohol or flavored drinks (e.g. fruit juice), breath sprays containing alcohol, medicinal juices and drops can distort the result.
- Belching and vomiting can distort the result.
- Drafts and climatic conditions such as wind, moisture and dust can distort the result when performing the measurements with the funnel.
- If the device is used for maritime application, a distance of 5 m to the compass must be complied with.

4.4 Preparations for use

1. Ensure that the rechargeable battery is charged or that the batteries (for battery-operated devices) have been inserted.
2. Turn on the device by pressing the OK button.
 ⇒ A self-test is carried out.
 When the device is ready for measuring: Wait for the display to indicate that the device is ready to measure (see the fold-out page, Figure 2: **A**).

4.4.1 Before sampling

- The subject must be breathing calmly and steadily.
- The subject should be able to provide the minimum volume of breath (>1.2 L) required for sampling.

4.5 Changing the measuring mode

By default, the device is set to measuring using the mouthpiece. The measuring mode must be changed before measuring using the funnel.

ℹ The active measuring mode is indicated by a corresponding symbol (device with mouthpiece or funnel) on the display.

To change the measuring mode:

1. Press the ^ button to access the context menu.
2. Choose ***Funnel***.
3. Press OK to confirm.

✓ The Funnel measuring mode is now set.

ℹ The measuring mode remains set until it is changed in the menu or the device is restarted. Restarting returns the device to the mouthpiece measuring mode.

4.6 Breath sample

4.6.1 Taking an automatic breath sample

ℹ Use a new mouthpiece for each measurement.
Do not expose the sensor to high alcohol concentrations, e.g. by rinsing the mouth with alcoholic drinks with a high alcohol proof shortly before the measurement. Doing so shortens the life span of the sensor.

ℹ Before taking the breath sample, ensure that the correct measuring mode has been selected.

1. Remove the mouthpiece (see "Order list", page 54) from its packaging. For hygiene reasons, keep the mouth contact area covered by the packaging to protect it until the mouthpiece is firmly attached to the device.
2. Position the mouthpiece at an angle with the serrated marking (contact surface) facing upward and push it up and into the mouthpiece receptacle (see Figure 1 on the fold-out page).
 ⇒ The mouthpiece audibly clicks into place.

 The mouthpiece can be used at a left or right angle.
3. Fully remove the packaging.
4. Blow evenly and without interruption into the mouthpiece.
 ⇒ A circular progress bar indicating progress as a percentage value is displayed (see the fold-out page, Figure 2: **B**, **C**).

4.6.2 Reading the measurement result for automatic breath samples

The following settings depend on the configuration.
The pre-configured alcohol threshold values determine when the green, yellow or red LED illuminates.

Three types of results are possible.

- No alcohol detected (see the fold-out page, Figure 2: **D**). In addition, the green LED lights up.

Or

- Low level of alcohol detected (see the fold-out page, Figure 2: **E**). In addition, the yellow LED lights up.

Or

- Alcohol detected (see the fold-out page, Figure 2: **F**). In addition, the red LED lights up.

Press the OK button to continue with the next measurement.

Once the recovery time of the sensor (see "Technical data", page 51) has elapsed, the display indicates that the device is ready for measurement.

4.6.3 Taking a manual breath sample (optional)

If the subject cannot provide the required minimum volume of breath, manual sampling can be activated.

To take a manual breath sample:

1. Turn on the mode under breath sample options in the PIN-1 menu.
2. Follow steps one to four of the process for taking an automatic breath sample (see "Taking an automatic breath sample", page 37).
3. If the progress circle is displayed after starting the breath sample, tap the OK button.
 ⇒ The following note appears in the results display: ***MANUAL***
4. Press the OK button to continue with the next measurement.

4.6.4 Taking a passive breath sample

During passive sampling, the sampling is manually initiated by the device user. Possible applications for this include testing the ambient air or a person's breathing air for the presence of alcohol.

To take a passive breath sample:

1. Hold the device, which is ready to measure, into the air to be analyzed **without** a mouthpiece.
2. Tap the OK button to initiate the measurement.
 ⇒ If alcohol is detected, the message ***ALCOHOL*** is displayed and the red LED lights up (see Figure 2 on the fold-out page: **H**).
 ⇒ If no alcohol is detected, the message ***NO ALCOHOL*** is displayed and the green LED lights up (see Figure 2 on the fold-out page: **G**).
 ⇒ In the results display, ***PASSIVE*** is also displayed.
3. Press the OK button to continue with the next measurement.

4.7 Breath sampling with the funnel

4.7.1 Taking a breath sample

⚠ CAUTION
Hygiene

- ► Do not touch the upper part of the funnel.
- ► Always maintain a distance of 2-3 cm from the object or subject to be tested.
- ► Do not touch excretions, surfaces or a subject's skin with the funnel.
- ► Avoid direct contact with substances to be tested.
- ► Replace contaminated funnels.

⚠ CAUTION
Incorrect measurement results
If multiple samples are taken over a short period of time, there is a risk of condensation building up inside the funnel and entering the device.

- ► Replace the funnel if condensate or droplets are visible.

ℹ Before taking the breath sample, ensure that the correct measuring mode has been selected.

1. Remove the funnel from its packaging.
2. Hold the funnel using the grip areas on the sides and attach it firmly to the top of the device (it can only point in one direction).

3. Ask the subject to blow evenly and without interruption into the funnel without touching it.
 ⇒ A continuous tone signals a sufficient flow of air.
 ⇒ While the sample is being taken, an animated circular progress bar that indicates progress as a percentage value is displayed (see Figure 2 on the fold-out page: **B**, **C**).
 ⇒ If sampling was successful, the tone ends with a clear, audible click.
 ⇒ The display shows a spinning progress bar during an ongoing analysis process.
 ⇒ The result is displayed after a few seconds and stored by the device.

4.7.2 Reading the result

- If alcohol is detected, the message ***ALCOHOL*** is displayed and the red LED lights up.
- If no alcohol is detected, the message ***NO ALCOHOL*** is displayed and the green LED lights up.
- Press the OK button to continue with the next measurement.

4.8 Turning off the device

- Hold down the OK button for about 3 seconds.

The device will switch off automatically approx. 4 minutes after the last use.

4.9 Removing the mouthpiece

- Push the mouthpiece on the mounting spacer out of the mouthpiece receptacle in an upward motion and dispose of it (see fold-out page, Figure A).

4.10 Removing the funnel

- Hold the funnel using the grip areas on its sides, pull it out of the device in an upward motion and dispose of it.

4.11 Menu items

The menu (see "Menu and context menu", page 35) contains the following items:

4.11.1 Test results

Displays all stored results with the measuring mode, measuring unit, date, time and sample number, including whether the sample was taken manually or passively.

4.11.2 About screen

Displays the part number and serial number of the device and the EC sensor.

4.11.3 Next service

Displays the time until readjustment or service is due in days.

4.11.4 Bluetooth®

NOTICE
Risk of other devices malfunctioning.
This device includes the Bluetooth® feature. Bluetooth® must not be used on aircraft.

► Do not select Bluetooth® on aircraft.

The device is ready to print via Bluetooth® and to connect to an app via Bluetooth®.
An active Bluetooth® connection is indicated by the Bluetooth® symbol on the display.

Turns the Bluetooth® feature on or off. Displays coupled and last used printers. Searches for available printers and removes printers.

Bluetooth on/off:

- Choose Bluetooth® on or off with the or button.
- Press OK to confirm.

Recent devices:

- Choose the desired printer using the and buttons.
- Press to confirm.

⇒ The Bluetooth® connection to the printer is established. The printer appears in the menu.

Scan devices:

1. Press the button.
 ⇒ Available printers are searched for and displayed.

To couple and connect a printer for the first time, see the Dräger Mobile Printer BT instructions for use.

Remove devices:

1. Choose the desired printer using the and buttons.
2. Press to confirm.
3. Press the button to confirm the query.
 ⇒ The printer is removed.

4.11.5 Date and time

Sets the format, date and time.

- Choose the date or time format using the or button.
- Press to confirm.
- Change the digits using the or buttons.
- Press to confirm.

4.11.6 Brightness

Adjusts the brightness of the display screen.

- Adjust the brightness using the and buttons.
- Press to confirm.

4.11.7 Sounds

Turns the sounds on or off.

- Turn the sounds on or off using the and buttons.
- Press to confirm.

4.11.8 Language

Sets the language.

- Choose the language using the and buttons.
- Press OK to confirm.

4.11.9 Measuring unit

Sets the measuring unit.

The default unit setting for the measured value of breath alcohol is mg/L (= milligrams of alcohol per liter of breath). The device can be set to all common measuring units.

- Choose the measuring unit using the and buttons.
- Press OK to confirm.

4.11.10 Automatic turning off

Sets the time until automatic shutdown.

- Set the time using the and buttons.
- Press OK to confirm.

4.11.11 Battery type of the AA batteries (only for AA variant)

Defines the battery type of the AA batteries that are used in the device.

- Choose the battery type using the and buttons.
- Press OK to confirm.

4.11.12 Software version

Displays the software version number and the associated part number.

4.11.13 Sign in

Access to other menu items.

- Enter the four-digit PIN-1.
- Confirm the entry twice with the OK button.

ℹ Access to other menu items is protected by entering a PIN.
The default PIN is: 2 0 2 0.
Dräger recommends changing the PIN on first use.

4.11.14 Turn off

Turns off the device.

- Tap the OK button.

⇒ The device turns off.

4.12 GPS (Global Positioning System, optional)

NOTICE
Risk of other devices malfunctioning!
If GPS is installed, the device must not be used on aircrafts.

► Do not use the device when on aircraft.

The unit can be equipped at the factory with GPS (Global Positioning System) for the satellite-based determination of the device's position.

If the GPS option has been installed, the determined position data is stored together with the measurement results and output when the device is read out.

4.13 Printing Results

To print the measurement results, the device can be connected to the Dräger Mobile Printer BT (see "Order list", page 54) via Bluetooth®.

ℹ Observe the Dräger Mobile Printer BT instructions for use.

Alternatively, the measurement results can also be printed using the Dräger Mobile Printer.

ℹ Observe the Dräger Mobile Printer instructions for use.

- In the menu, select ***Login***.
- Enter PIN-1.
- Press OK to confirm.

- In the PIN-protected menu, select ***Printer***.
- Under ***Printing mode***, select the Mobile Printer or Mobile Printer BT.

4.14 Exporting data

The data stored by the device can be exported and then analyzed in Microsoft® Excel, for example. (This function is available from firmware version 2.0.0).

1. Connect the USB-C cable to the USB-C port of the device and the USB port of a PC or laptop.
 ⇒ The device turns on.
2. Log in to the PIN-1 menu in the device menu and choose ***Transfer mode***.
 ⇒ "***Transfer mode***" is shown on the device display.
3. Open Windows Explorer.
 ⇒ The device is now displayed as a USB drive.
4. Open the path “USB drive\var\scenario_logs\BreathTest”.
5. Select the CSV file.
6. Select Copy.
7. Select the desired drive and directory as the storage location.
8. Select Save.
 ⇒ The CSV file is saved.
9. Turn off the device by pressing the OK button to exit ***Transfer mode***.

4.15 App connection

The device data can be transferred to the Dräger add app automatically via Bluetooth.

After initial pairing, the Dräger add app connects to the device automatically and downloads the data. The data are displayed there as a list, in a ring graphic and also on a map if GPS is activated. Other data can be added in addition to the measurement results.

For initial pairing:

1. Install the Dräger add app by following the instructions on the screen.
2. Choose the device in the app.
3. Choose Bluetooth in the device menu.
4. Select Bluetooth and choose Pair with app.
5. Confirm the connection in the app.
6. Compare the code on **both** devices and confirm.

5 Troubleshooting

An error message is displayed in the event of a device malfunction. The number below the message is for service purposes. If the error persists even after turning the device off and on multiple times, contact authorized personnel, such as DrägerService.

Error/Message	Cause	Remedy
Device does not turn on.	(Rechargeable) batteries are empty.	Charge the rechargeable battery or replace the batteries.
Insufficient volume of breath/ invalid sample.	Subject is not blowing for long enough or intermittently.	Blow more strongly and without interruption into the mouthpiece/funnel. If necessary, take a manual or passive breath sample.
Device out of temperature range.	Device is too warm or too cold.	Allow the device to cool down or warm up.

Error/Message	Cause	Remedy
↑ ↑ ↑ (Measuring range exceeded).	Measuring range exceeded. Residual mouth alcohol can distort the measured values.	A wait time of at least 15 minutes after alcohol was last consumed must have elapsed. To repeat, press the OK button.
Insufficient volume of breath. [Only for breath samples using the mouthpiece]	Volume of breath sample insufficient.	Repeat the breath sample after pressing the OK button. If necessary, take a manual or passive breath sample.
Sampling interrupted.	Irregular exhalation, abrupt cessation or inhalation.	Repeat the breath sample after pressing the OK button. If necessary, take a manual or passive breath sample.
Error number.	Device error.	Turn the device off and then back on. If the error persists, contact DrägerService or your dealer.
The red LED flashes three times in close succession and the device turns off.	(Rechargeable) batteries are empty. Device is too warm or too cold.	Charge the rechargeable battery or replace the batteries. Allow the device to cool down or warm up.

6 Maintenance

The device is delivered with default adjustment settings. The date of the last adjustment can be recalled in the menu.

Dräger recommends that the device be adjusted once a year. Adjustment should only be carried out by authorized persons, e.g., DrägerService.

Calibration certificates:
https://www.draeger.com/productioncertificates

ℹ Comply with the national guidelines governing adjustment.

6.1 Charging the rechargeable battery

⚠ WARNING
Risk of burns and electric shock!
The use of power supply units that have not been approved by Dräger can cause burns and electric shocks when charging the rechargeable battery.

► Use only the power supply unit specified in the order list.

- Connect the USB-C cable to the USB-C port of the device and a USB charger.

⇒ The rechargeable battery will charge.

6.2 Replacing the Batteries

⚠ WARNING
Risk of explosion
Risk of explosion if the incorrect battery type is used.

► Only use the battery types specified by Dräger.

1. Turn off the device.
2. Slide off the battery cover.
3. Replace the batteries (3 x 1.5 V alkaline battery [LR06 Mignon AA]), ensuring that they are inserted the right way round.
4. Slide the battery cover back into place.

ℹ The time must be reset if the batteries are not replaced in less than 5 minutes. Measurements cannot be performed if the time has not been set.

6.3 Cleaning

Clean the device using an alcohol-based cleaning agent.

To clean the device:

1. Switch off the device.
2. Remove the mouthpiece from the device.
3. Wipe the device down with a disposable cloth moistened with detergent. Do not apply the cleaning agent directly to the device.

 ℹ Ensure that no liquid gets into the metal piece of the mouthpiece receptacle on the top of the device.

4. Check that there is no residual alcohol stuck to the device via passive sampling and that the result ***NO ALCOHOL*** is displayed.
5. If ***ALCOHOL*** is still displayed, repeat passive sampling after some time (10–15 minutes) until the result ***NO ALCOHOL*** is displayed.
6. If ***ALCOHOL*** is still displayed, store the device for up to 24 hours to allow the detergent to evaporate. Repeat passive sampling after some time.

For information on disinfecting Alcotest devices, see document 9300698 at www.draeger.com/IFU.

7 Storage

- Do not permanently store the device in its delivery packaging.

⚠ CAUTION
Risk of damaging the device!
Leaking battery acid can damage the device.

► Remove the batteries if the device is not used for a prolonged period of time.

Maintain the charge level of the rechargeable batteries.
For nickel-metal hydride batteries (charge outside the device):
To ensure operational readiness of the device, even if it is not used, fully charge the rechargeable battery at least every three months.
For lithium-ion batteries:
Charge to at least 50% every 6 months if not used for a prolonged period of time to avoid shortening the life span.

8 Disposal

This product must not be disposed of as household waste. This is indicated by the adjacent icon.
You can return this product to Dräger free of charge. For information please contact the national marketing organisations and Dräger.

Batteries must not be disposed of as household waste. They are therefore marked with the adjacent symbol. Dispose of batteries and rechargeable batteries at battery collection centers as specified by the applicable regulations.

9 Technical data

Measurement principle	Electrochemical sensor
Measuring range	
Breath alcohol concentration	0.00 to 3.00 mg/L (mass of ethanol per minimum volume of breath at 34 °C and 1013 hPa)
Blood alcohol concentration	0.00 to 6.00‰ (mass of ethanol per blood volume or mass of blood at 20 °C and 1013 hPa)
A conversion factor of 2100 is used to calculate the blood alcohol concentration using the breath alcohol concentration data (factory setting, country-specific conversion factors can be configured).	
Measuring units	

mg/L, µg/100 mL, µg%, µg/L, g/210 L, g/L, ‰, %, mg/100 mL, mg%, mg/mL	
Ambient conditions	
During operation	–10 to +55 °C 10 to 100% r. h. (non-condensing) 600 to 1300 hPa
In storage	–20 to +60 °C 15 to 75% r. h. (10 to 100% r. h. <48 h)
Data storage	
	Up to 10,000 tests (depending on the configuration) Warning that data storage is almost full at 90% Preset: FIFO (oldest test is overwritten) Optional: Device locks at 100%
Sampling	
Minimum volume of breath	>1.2 L (default)
Breath flow	>10 L/min (default)
Breath duration	>1.6 s (default)
Measurement accuracy	
Standard deviation of repeatability using standard alcohol	
up to 0.50 mg/L ≥ 0.50 mg/L up to 1.00‰ ≥ 1.00‰	0.008 mg/L 1.7% of the measured value 0.017‰ or 1.7% of the measured value Depending on which value is higher
Sensitivity drift	Typically 0.4% of the measured value per month
Wait time until next measurement can be taken following a previous measurement of:	

0.00 mg/L 2 s	0.00‰ 2 s
0.25 mg/L 20 s	0.50‰ 20 s
0.50 mg/L 40 s	1.00‰ 40 s
1.00 mg/L 80 s	2.00‰ 80 s
>1.50 mg/L 120 s	>3.00‰ 120 s
Adjustment	Recommendation: every 12 months
Dimensions	61 mm × 172 mm × 36 mm (W × H × D)
Weight	
AA variant	285 g
Li-Ion variant	275 g
Power supply	
Lithium-ion battery	Permanently installed in the device
Lithium primary cell	3 x 1.5 V (Mignon AA), not rechargeable
[1)]Nickel-metal hydride battery	3 x 1.2 V (Mignon AA), not rechargeable via device
[1)]Alkaline batteries	3 x 1.5 V alkaline (LR06 Mignon AA)

Lithium-ion battery for up to 5000 measurements,
batteries for up to 1500 measurements using the mouthpiece/funnel (depending on use and ambient conditions)

Interface	USB-C
LED risk group	2
Protection class	IP 54
CE marking	Electromagnetic compatibility
	Low Voltage Directive

The device meets the requirements of EN 15964. It can be configured in accordance with the standard. Display of measured values if configured in accordance with EN 15964:

Measured value	**Display**
0.0 mg/L to 0.03 mg/L	0.0 mg/L
> 0.03 mg/L	Measured value

1) At low temperatures, the number of tests may be limited (EN 15964).

10 Order list

Designation and description	Order no.
Mouthpiece with exhalation non-return valve (100 pieces)	6811055
Mouthpiece with exhalation non-return valve (250 pieces)	6811060
Mouthpiece with exhalation non-return valve (1000 pieces)	6811065
Mouthpiece without exhalation non-return valve (100 pieces)	6810690
Mouthpiece without exhalation non-return valve (250 pieces)	6810825
Mouthpiece without exhalation non-return valve (1000 pieces)	6810830
Funnel (10 pieces)	8327718
AA lithium battery	3700481
AA alkaline battery	1335804
AA NiMH battery	1890092
Wrist strap	3700416
USB-C cable	3700418
Vehicle USB charging adapter	3730475
Power supply unit	3701090
Dräger Mobile Printer BT	3700421

Manufacturer
Dräger Safety AG & Co. KGaA
Revalstraße 1
D-23560 Lübeck
Germany
+49 451 8 82-0

Importer (UK):
Draeger Safety UK Ltd
Ullswater Close
Blyth Riverside Business Park
Blyth, Northumberland NE24 4RG
United Kingdom

9300088 – 4754.591 me
© Dräger Safety AG & Co. KGaA
Edition: 08 – 2025-01 (Edition: 1 – 2020-10)
Subject to alterations

www.draeger.com